



*MINI*

HYPERSPHERE

## **IMPORTANT SAFETY INSTRUCTIONS**

When using an electrical appliance, basic precautions should always be followed, including the following:

**READ ALL INSTRUCTIONS BEFORE USING THE HYPERSPHERE MINI**

### **DANGER**

To reduce the risk of electric shock:

- Do not reach for an appliance that has fallen into water. Unplug immediately.
- Do not use while bathing or in a shower.
- Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

### **WARNING**

To reduce the risk of burns, fire, electric shock, or injury to persons:

- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- Do not carry this appliance by supply cord or use cord as a handle.
- Keep the cord away from heated surfaces.
- Do not use outdoors.

- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not immerse unit in water. Keep liquids away from ventilation ports, buttons and charging port.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. After one hour use, allow unit to rest for 30 minutes before re-using.
- Unplug the unit after charging or prior to use.
- Do not use directly on a hard surface. Use only on a carpeted floor or a fitness mat.
- Persons weighing above 250 lbs/113 kgs may damage or break the unit during normal intended use.
- This is not a toy. For adult use only. Do not use if injured. Consult your doctor before using this product.
- Use product only according to Instructional Use materials available at [www.hyperice.com](http://www.hyperice.com).

## **CAUTION**

The battery used in this device may present a risk of fire or chemical burn if mistreated. Do no disassemble, heat above 60° Celsius, or incinerate.


Use only with any standard 2.0 USB connection.

## **SAVE THESE INSTRUCTIONS**

## HYPERSPHERE MINI

The Hypersphere Mini is a cutting-edge vibrating massage ball that uses pressure and vibration to relieve tension and keeps muscles loose, and promotes increased range of motion and flexibility.



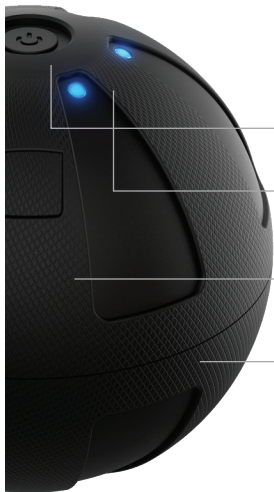

## SPECIFICATIONS

- 1** Compact Size – 3.9" (8.9cm)
- 2** Digital Circuitry Controls  
w/ 3 Speed Settings
- 3** USB Rechargeable Lithium Ion  
Battery w/ Power Indicator  
Lights (2+ Hr of Use Per Charge)
- 4** High-Intensity Vibrating  
Core Transfers Maximum  
Vibration to the Body
- 5** Textured Rubber  
Exterior Targets Vibration  
For More Localized Therapy

## USAGES

The Hypersphere Mini helps to loosen the body's soft tissue (muscles, tendons, ligaments, fascia). With consistent use, the Hypersphere Mini may relieve tension and keeps muscles loose, and promotes increased range of motion and flexibility. With its compact size, the 'Mini' is perfect for travel and on-the-go.



AT WORK



IN THE CAR



AT HOME



TRAVELING

## FEATURES

- Digital circuitry control
- 3 speeds (42Hz / 52Hz / 62Hz)
- High intensity vibrating core
- Patented design with textured silicone exterior
- Rechargeable lithium-ion battery (2+ hrs of use per charge)
- Lightweight, compact size (1 lb / 3.5")
- TSA approved for carry-on

## IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING THE HYPERSPHERE MINI

Dispose of used battery promptly. Keep away from children. Do not disassemble and do not dispose of in fire.

**Please do not use the Hypersphere Mini, or any vibration high-intensity exercise device without first obtaining approval from your doctor if any of the following apply:**

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any vibration device.

These contra indications do not mean that you are not able to use a vibration or exercise device but we do advise you to consult a doctor first.

Ongoing research is being done to the effects of vibration exercise for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above. Practical experience has shown that integration of vibration exercise into a treatment plan is advisory in a number of cases. This must be done on the advise of and in accompaniment with a doctor, specialist or physiotherapist.



## **OPERATING INSTRUCTIONS**

READ ALL INSTRUCTIONS BEFORE USING THE HYPERSPHERE MINI

### **CHARGING**

Fully charge for two hours before first use. To charge, connect the USB cable from the charging port to any standard USB charger (not included). LED lights illuminate to show battery level and indicate active charging. LED lights will remain illuminated when fully charged while adapter is plugged in. Unit may be recharged at any time and at any battery level. It is not recommended to fully deplete the battery to the red LED level. Recharge the unit monthly to best prolong battery life. The average run-time is 2 hours at level three, and longer at lower levels.

### **OPERATING**

To turn on, press down on the power button for three seconds until LED lights illuminate. Press again once for level one, twice for level two and three times for level three. The blue LED lights will indicate the corresponding power level. Do not use directly on a hard surface. Use only on a carpeted floor or a fitness mat. To turn off, press down on the power button for three seconds. The unit is fully turned off when all the LED lights are off. The unit will turn off after 10 minutes of use. To continue use, turn unit on again and set to desired level.

## OPERATING INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING THE HYPERSPHERE MINI


- 1 Power Button
- 2 Battery level indicator LED (green, red)
- 3 Frequency level indicator LED (blue)
- 4 USB Charging Port
- 5 USB Charging Cable (not shown)

### SPECIFICATIONS

- Rated Charging Voltage: 5V DC- 1A
- Rechargeable Lithium Ion battery 2400mAh
- Frequency: Low 42Hz, Medium 52Hz, High 62Hz
- Motor RPM: Low 2500, Medium 3100, High 3700
- Weight: .75 lbs / .34 kg
- Designed in USA. Assembled in China

### MAINTENANCE, CLEANING & STORAGE

- Use a towel or soft cloth to wipe any moisture from unit.
- Keep unit away from heat sources, liquids, and children.

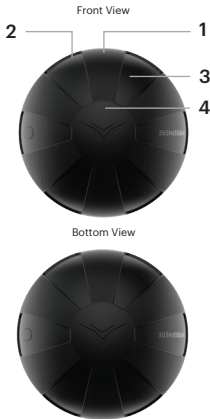

## 1 YEAR WARRANTY

**If within one year from the date of purchase, this product fails due to a defect in material or workmanship, Hyperice, Inc. will repair or replace the product, or necessary components, free of charge.**

This warranty excludes:

- (a) damage caused by accident, abuse, mishandling, or transport;
- (b) units subjected to unauthorized repair;
- (c) units not used in accordance with Hyperice care instructions;
- (d) damage exceeding the cost of the product;
- (e) deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) failure to provide the dated proof of purchase.

Some states and countries do not allow a limitation of damages, so the foregoing limitation may not apply to you. This warranty guarantees specific legal rights, and but other rights may vary from country to country and from state to state. This warranty service is available through the operating policies and procedures of Hyperice Inc. and all Hyperice international distributors. If you are a client, please contact the appropriate Hyperice business entity, dealer or retailer, from whom you directly purchased the product for warranty and return authorization procedures. This warranty is initiated and executed by Hyperice Inc.

## OTHER PRODUCTS

### Hypervolt™

#### Vibration Massage Device

Featuring Quiet Glide™ Technology

The Hypervolt is a cordless state-of-the-art vibration massage device that helps relax sore and stiff muscles to improve mobility.



### Vyper™ 2.0

#### High-Intensity Vibrating Fitness Roller

The Vyper 2.0 is the next generation of the world's most powerful vibrating fitness roller. Now stronger, sleeker, and even more durable.



**Available Colors:** Black, Green, Grey Camo, Red Camo

## OTHER PRODUCTS

### Hypersphere™

#### High-Intensity Vibrating Fitness Ball

Powerful yet small, the Hypersphere targets your tightest areas to unlock the body and help you move better.



**Available Colors:** Green, Black

### Venom™

#### Wearable Vibration + Heat

The Venom is a cutting edge, digitally connected wearable device that combines heat and vibration to warmup, loosen and relax muscles.



**Available Devices:** Back, Right Shoulder,  
Left Shoulder, Leg

## OTHER PRODUCTS

### Hyperice® ICT

Ice Compression Technology  
w/ Patented Air Release Valve

The product that started it all. Designed in collaboration with the world's best athletes, this recovery staple paved the way for the brand's future innovations.



**Available Devices:** Right Shoulder,  
Left Shoulder, Knee, Back, Utility



This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.



525 Technology Drive, Suite 100

Irvine, California 92618

(949) 565-4994

**hyperice.com**

© 2021 Hyperice, Inc. All Rights Reserved.

170-00001 Rev 1